

# Sure-Lok®

## Needle-free IV Connector



Part of Amsino's AMSafe® product line, the Sure-Lok® is designed with transparent housing, tight septum seal and smooth swab able surface to help reduce bacterial contamination and support infection control.<sup>1</sup>

The Sure-Lok® has a high flow rate and a proven lower resistance to luer connections allowing ease of connection the ability to use the device easily and effectively.<sup>2</sup>

### Features – Benefits


- **Clear housing and septum** – Provides a visual confirmation of fluid flush
- **Improved flow rate** – Fluid delivery at a rate >200 mL/minute
- **Lower resistance to luer connections** – Designed for easy connection and disconnection
- **Flat top with tight septum seal** – Helps to support effective disinfection
- **Metal free** – MRI compatible

### Standard IV Extension Sets

REF	LENGTH	DESCRIPTION	APPROX. PRIMING VOLUME	QTY/UNIT
<b>AE3108</b> ©	8 in. (20 cm)	Removable Sure-Lok Needle-free Connector, Pinch Clamp	1.2 mL	100/Case
				
<b>AE3108B</b> ©	8 in. (20 cm)	Bonded Sure-Lok Needle-free Connector, Pinch Clamp	1.2 mL	100/Case
				

### Sure-Lok Needle-free Connector

REF	DESCRIPTION	APPROX. PRIMING VOLUME	QTY/UNIT
<b>AY0200</b> ©	Sure-Lok Needle-free Connector	0.33 mL	100/Case

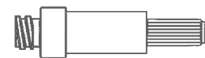

View our entire line of products at [www.amsino.com](http://www.amsino.com)



#### Amsino Medical Group Co. Ltd.

708 Corporate Center Drive, Pomona, CA 91768 USA  
+1 800.632.6746 • fax: +1 909.626.3888 • [info@amsino.com](mailto:info@amsino.com)

#### SOURCES:

1. Jarvis W, MD. Choosing the Best Design for Intravenous Needleless Connectors to Prevent Bloodstream Infections. Infection Control Today, August 2010 <https://www.infectioncontroltoday.com/infusion-vascular-access/choosing-best-design-intravenous-needleless-connectors-prevent-bloodstream>
2. Amsino International, Inc. has proprietary data on file.

Product(s) described may not be licensed or available for sale in all countries. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. AMSafe and Sure-Lok are registered trademarks of Medical Group Co. Ltd. © Copyright 2024, Amsino Medical Group Co. Ltd. All rights reserved.